

## 4 Way Flow Selector and Frame

### SPECIFICATION SHEET

With BS Schrader Input Probe on 2m hose.


#### Features

MEC Medical Ltd Flow Selectors are designed and built to the highest standards with clear display of the flow rate.

- Strong aluminium frame.
- Extending hooks for hanging.
- 2m supply hose to connect to 4 bar oxygen source.
- Folding carrying handle.
- Four therapy outlets, each with selectable flows, ½, 1, 2, 3, 4, 6, 8, 10 and 15 litres per minute.
- Plastic feet to base.

#### Specification

Weight including hose:	3.4Kg
Length:	46cm
Height:	18.5cm
Height with handle raised:	24cm
Depth:	15.5cm

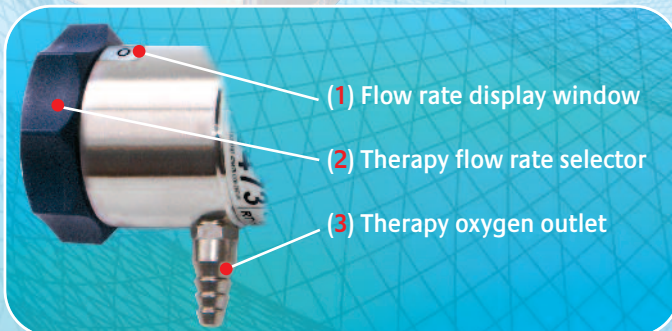

Oxygen therapy is easily chosen by rotating the "Therapy flow rate selector" (2). The flow rate is indicated in the display window (1). Choose from ½, 1, 2, 3, 4, 6, 8, 10 and 15 litres per minute. Oxygen at the chosen flow is discharged from the Therapy outlet (3). The Therapy outlet is designed so that plastic oxygen tubing can be securely held in place on the fir-tree style connector.

**Note:** The flowmeters are designed so that the flow rate display window can be easily seen, at the same time the therapy outlet can be accessed so that equipment may be attached.



Supplied by